①

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

ⒼⒷ

This is a special support for the so called tennis elbow. Pain in the elbow will occur through overstraining. Application of this brace 2" below the elbow joint of the forearm will reduce the tension in this area.

#### Instructions for use:

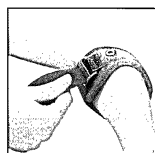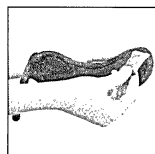Thread the velcro strap into plastic ring and then slip the strap onto your forearm. To position correctly, please take care of the marking, printed on the bandage. The narrow part of pelotte should always show in the direction of elbow joint. The inserted pelottes are to exert pressure on the muscle elevation from the back of the hand to the elbow. You can feel this muscle by squeezing your fist. The pressure of the BORT Epicondylitis Strap should only be felt on this muscle when you squeeze your fist or grasp something with your hand. If the hand is re-laxed, the strap should fit without any pressure being felt from the pelottes. Close the velcro fastening to make the strap fit tightly without hindering the blood circulation. Cut off the overlapping velcro strap.

#### Indications:

For pains in the forearm muscles near the elbow caused by sport injuries, excessive strain and „tennis-elbow“.

#### Important Instructions:

No reutilization – this medical device is meant for the care of one patient only.

Wear bandage only when under strain, not during rest or sleep. Avoid hindering blood circulation. In case of numbness loosen bandage and remove it if necessary. If complaints continue, please contact the doctor.

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